


panmed | **u.s.**
PAN | MEDICAL

Endovascular

Remodelling made easy



Your partner in healthcare™

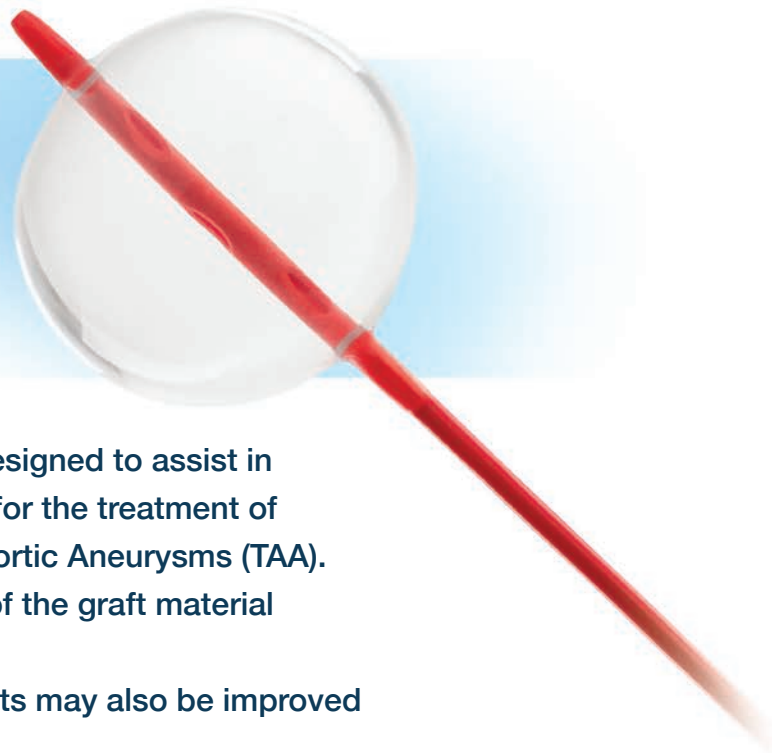

Constellation

Stent graft balloon catheter

The Constellation Stent Graft Balloon Catheter is designed to assist in the expansion of a self-expanding stent graft used for the treatment of Abdominal Aortic Aneurysms (AAA) and Thoracic Aortic Aneurysms (TAA). Ballooning of a stent graft may improve modelling of the graft material and fixation of the stent graft to the vessel wall. Sub-optimal expansion of self-expanding stent grafts may also be improved by inflating the balloon at the site of the stent.

Features

- Strong, conformable balloon.
- Broad range of diameters, ideal for AAA procedures.
- Low profile.
- Radiopaque marker bands at each end of balloon.
- Rapid inflation and deflation with multi-lumen placement catheter.
- Compliant polyurethane balloon.
- Extension with 3-way stopcock.

Inflation Volumes

Diameter (mm)	ml (cc)
10	3
20	9
30	19
40	41
46	60

Product Codes

Product Code	Catheter Length (cm)	Balloon Diameter (mm)
CT-RMB	100	10 - 46